

## Set de verificación de calibración de cloruro Checker® HC

HI753-11



## Descripción

El set de verificación de calibración de cloruro Checker® HC **HI753-11** proporciona una solución simple para validar el cloruro checker® HC **HI753** que tiene un rango de 0,0 a 20,0 ppm de cloruro. Este conjunto de estándares de alta calidad se fabrica en nuestras instalaciones de última generación y se suministra con un certificado de análisis. El certificado de análisis proporciona el número de lote, los valores de referencia y la fecha de vencimiento para la trazabilidad al certificar el Checker®HC.

## Especificaciones

No Especifica

## Accesorios

No Especifica

## Cómo pedir

No Especifica

## Ventajas

El **HI753-11** es un conjunto de estándares de calibración de alta calidad que permite a los usuarios validar el cloruro Checker<sup>®</sup> HC **HI753**. El **HI753-11** incluye una cubeta sellada para cero ppm de cloruro y una cubeta teñida a una muestra de cloruro equivalente a 10,0 +/- 1,0 ppm a 25 °C.

**-Reutilizable con estabilidad de larga vida útil.**

**-Se suministra con certificado de análisis:**

- Código de producto
- Numero de lote
- Fecha de caducidad



## CERTIFICATE OF ANALYSIS

Product Name: CALIBRATION CHECKING SET FOR HI 711  
 Product Code: HI 711-11  
 Lot number: SC0372/14  
 Expiration date: October 2017

| Standard cuvette | Lot number | Standard Value @ 25°C (ppm) | Passed |
|------------------|------------|-----------------------------|--------|
| A, HI 711-11     | 2038       | 1.00 ± 0.05                 | ✓      |
| B, HI 711-11     | 2039       | 1.00 ± 0.05                 | ✓      |

Hanna Instruments certifies these standards meet the stated tolerance limits as indicated. These should be used for quality control to validate instrument function on the specified Hanna instrument only. Store standards in the storage box upright at room temperature. Do not use standard cuvettes.


Instructions for use:

1. Thoroughly wipe standard cuvette with HI731318 lint free cloth and only handle using top of cuvette. The cuvette wall must be dry and free of smudges.
2. At the C.1 prompt place A, HI 711-11 standard cuvette into the holder of the Checker and press operations button.
3. At the C.2 prompt place B, HI 711-11 standard cuvette into the holder and press the operations button.

Acceptable results must be between 0.95 and 1.05 ppm @ 25°C.

File number: CERT711-11\_SC0372/14

QA manager: Bogdan Munteanu

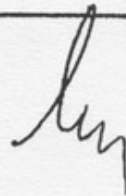

## Video

No Especifica